
RENPHO

User Manual

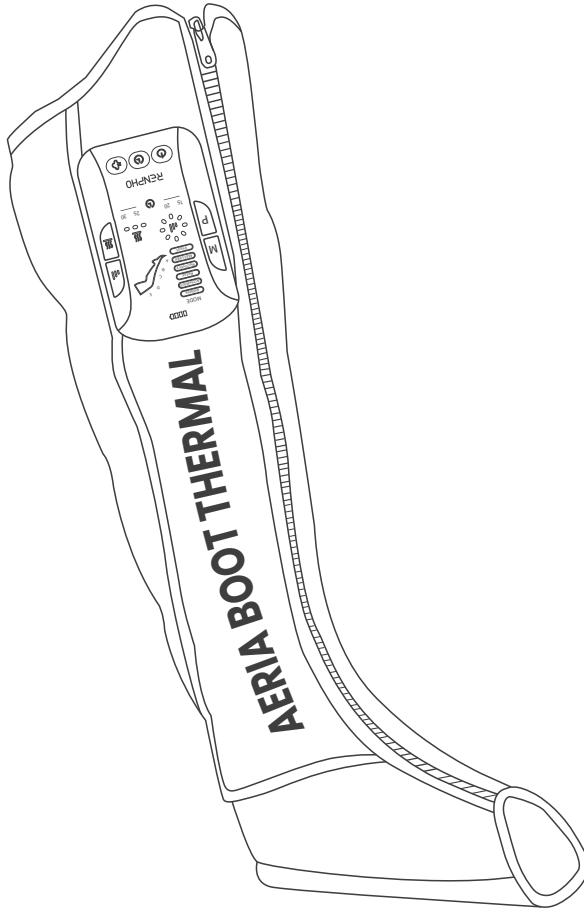AERIA BOOT THERMAL

Model: R-WML01

CONTENTS

• ENGLISH 01 ~ 10

• DEUTSCH 11 ~ 20

• FRANÇAIS 21 ~ 30

• ITALIANO 31 ~ 40

• ESPAÑOL 41 ~ 50

Thanks for purchasing the RENPHO Leg Massager. To ensure safety and avoid any accidental damage or injury, please read this manual carefully and use the device only as indicated. Non-compliance may cause injuries, and/or severe property damage.

IMPORTANT SAFETY INSTRUCTIONS

When using the device, basic precautions should always be followed, including the following:

READ ALL INSTRUCTIONS CAREFULLY BEFORE USING.

1. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the appliance.
3. Cleaning and user maintenance shall not be made by children.
4. The device is intended for household use only and not intended for medical use. It should not be used in treatment or alleviation of disease or abnormal physical state.

DANGER - To reduce the risk of electric shock and fire:

1. Never use pins or other metallic fasteners with this device.
2. The device is only to be used with the power supply unit provided with the device.
3. Do not leave the device unattended when it is working.
4. Always unplug the power adapter from the electrical outlet after each use.
5. Do not plug or unplug the power adapter or air hose with wet hands.
6. Keep the cord and device away from heated surfaces.

WARNING - To reduce the risk of fire, electric shock, or injury to persons, please follow the instructions below:

1. If you fall into any of the following user groups, please consult a doctor before use:
 - Pregnant or postpartum women
 - Children
 - Elderly (over 65 years old)
 - People with skin infections, sensitivity, redness, wound, and/or burns on the skin
 - People using a pacemaker or other implantable medical device
 - People with heart disease, malignant tumors, or other serious health conditions
 - People receiving medical treatment or feel discomfort during massages
 - People suffering from osteoporosis, diabetes, or thrombosis
 - People with high or abnormal blood pressure
 - People suffering from a sensory impairment that reduces the feeling of pain or heat (e.g. metabolic disorder)
 - People suffering from fever over 100.4°F
2. Stop using if you feel uncomfortable or unwell after use and consult a doctor.
3. **DO NOT** use spare parts that are not approved by RENPHO.
4. Non-professional repair personnel should not dismantle or replace the device without manufacturer's permission.
5. **DO NOT** use the device if the power cord or wrap surfaces are damaged.

6. **DO NOT** use or store the device in a place with high humidity such as bathroom or places with high temperature.
7. **DO NOT** carry the device by the power cord or the controller.
8. **DO NOT** use while sleeping, operating machinery, or after alcohol consumption.
9. **DO NOT** stand or walk while using the device.
10. **DO NOT** expose the controller, the wraps, or the power adapter to water as this may result in malfunction or electric shock.
11. **DO NOT** insert any foreign objects into the slit (this may cause fire and electric shock).
12. **DO NOT** use the sharp tools to scratch or puncture the wraps, as this may cause air leak.
13. **DO NOT** use wraps directly on swollen, inflamed areas, sensitive skin, or skin eruptions.
14. **DO NOT** use the device for more than 60 minutes at a time. Extensive use could lead to the device's excessive heating and shorter lifespan.
15. **DO NOT** operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury.
16. **DO NOT** use the device on pets or animals.
17. The device has a heated surface. Persons insensitive to heat must be careful when using the device.
18. **WARNING** - Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.

About the Battery

1. Do not expose the appliance or battery to excessive temperatures.
2. Be aware of the risk of terminals of the battery-operated appliance or battery being short-circuited by metal objects.
3. This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of.
4. Do not dismantle, open or shred battery.
5. Do not subject cells or batteries to mechanical shock.
6. In the event of a cell leaking, do not allow the liquid to come in contact with the skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical advice.
7. Secondary cells and batteries need to be charged before use. Always use the correct charger and refer to the manufacturer's instructions or product manual for proper charging instructions.
8. Do not leave a battery on prolonged charge when not in use.
9. After extended periods of storage, it may be necessary to charge and discharge the cells or batteries several times to obtain maximum performance.
10. **CAUTION** - The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not recharge, disassemble, heat above 100°C (212°F).

SAVE THESE INSTRUCTIONS

Battery Producer Labelling

Manufacturer: XIAMEN INVCY ELECTRONICS CO., LTD.



Address: 1F/3F-5F, BLDG 5, No.943-5, Tonglong 2nd Rd, Torch Hi-tech Industrial Zone (Xiang'an)Xiamen 361000 Fujian P.R. China

Website: <https://www.invcy.com/>

Email: 2162923595@qq.com

ABOUT THE PRODUCT


1. What's in the box

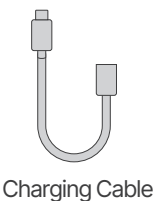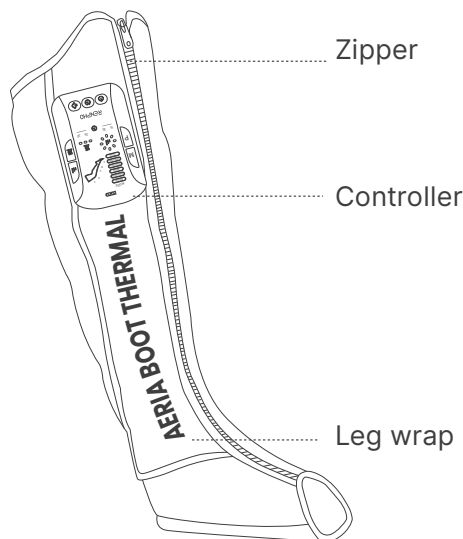Charging Cable

Buttons

- Power Button
- Timer Button
- Deflation Button
- M** Combine Mode Button
- P** Section Mode Button
- Intensity Button
- Heat Button

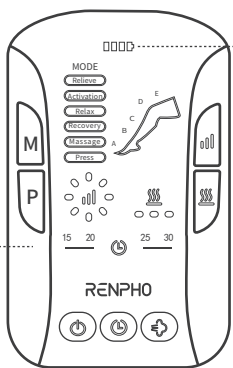LED Display

- Battery Indicator
- Combine Mode Indicator
- Section Mode Indicator
- Intensity Indicator
- Heat Indicator
- Timer Indicator

2. Specifications

Product Name: Leg Massager

Model: R-WML01

Input: 5V  2A

Battery Voltage: DC 7.4V

Battery Capacity: 4500mAh 33.3Wh

Auto-off Timer: 15/20/25/30 minutes

Charging Time: <360 min

Operating Time: about 160 min


Product Weight: 3.86 lb / 1.75 kg

Controller Size: 6.10×3.86×2.17 in / 15.5×9.8×5.5 cm

USING THE PRODUCT

1. Before Use


Charge the device

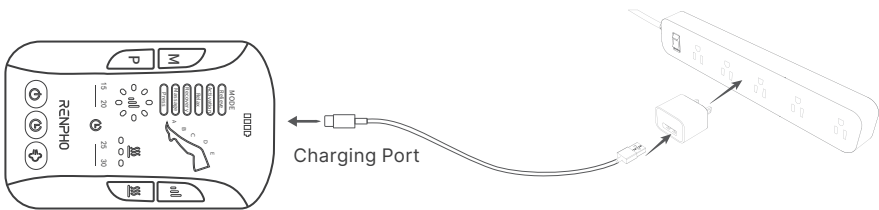
- a. Please charge the device using the provided USB-C charging cable and DC 5V 2A charger (not included) when the battery indicators turn off.
 - b. It takes about 6 hours to fully charge the device. When charging the battery, the battery indicator will flash and turn solid once the charging is complete. The massager runs on a single charge for about 160 minutes (varying by mode or intensity level selected).
- * For safety reasons, the device cannot be turned on while charging.

2. Wear the Leg Wrap

2.1 Lay the leg wrap flat.

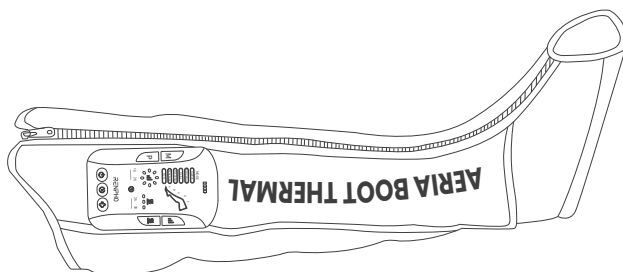2.2 Begin by unzipping the leg wrap. Insert your leg and re-zip to seal. The leg wrap is not size-adjustable and may provide a loose fit which is normal.

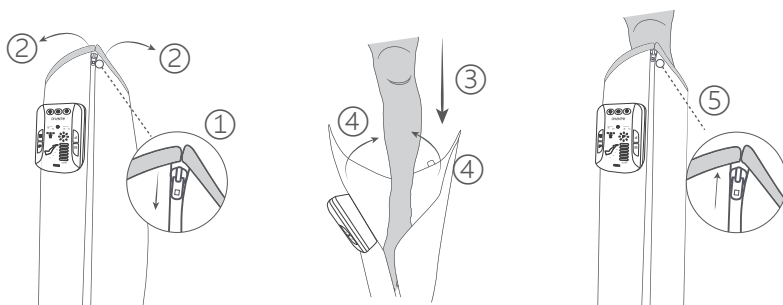2.3 Adjust your position and make sure your leg is fully extended.

3. Operation Instruction


1. **Power Button:** Press and hold to turn on / off the massager.
2. **Mode Button:** Press to toggle through 6 preset full-leg massage modes: Relieve-Activation-Relax-Recovery-Massage-Press.
3. **P Button:** Press to toggle 13 section modes of different combination, to target specific spot or multiple areas.
4. **Intensity Button:** Press to increase/decrease through 8 intensity levels (80mmHg to 220mmHg).
5. **Heat Button:** Press to toggle the heat level: Low, Medium, High and OFF.
6. **Timer Button:** Press to toggle massage duration: 15 min, 20 min, 25min and 30 min.
7. **Deflation Button:** Press to rapidly deflate the airbag for instant removal. This action also turns the massager off.

Note:

1. The massager has a built-in auto-off timer. It will automatically turn off the massager when the time is up. It is recommended to limit the usage time to 60 minutes per session.
2. If you experience any discomfort or if the leg wrap feels too tight, press and hold the Power button to turn off the massager immediately or press Deflation button for quick relief.
3. Memory function: The device automatically resumes the last selected massage mode upon power-up; Intensity, heat, and timer settings will revert to their defaults.

CLEANING AND MAINTENANCE

A. Storage

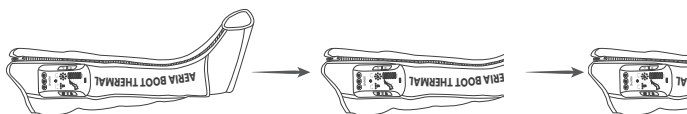Lay it flat and fold in half, then fold again.

1. Store the massager in a cool, dry place when not in use.
2. Avoid twisting the cord or folding the hoses when packing.
3. Do not place the massager in high-temperature or high-humidity locations. Avoid direct sunlight.
4. Avoid needles or sharp items that could puncture the airbags and hoses.
5. Do not place heavy items on the device.

B. Cleaning

1. Turn off the massager before cleaning.
2. Do not machine wash. Do not immerse it in any liquid to clean. Do not clean it with abrasive cleaners, brushes, gasoline, kerosene, glass polish, or paint thinner.
3. Wipe down the surface of the massager with a soft, slightly damp cloth instead.

C. Maintenance

1. Do not attempt to disassemble or repair the massager. The massager does not contain user-serviceable parts.
2. Keep the massager clean and away from moisture.
3. Cover up the massager when not in use to prevent dust.

TROUBLESHOOTING GUIDE

PROBLEM	CAUSE	SOLUTION
The device does not work.	The device is not turn ON correctly.	Press and hold the power button to turn on.
	The built-in auto-off timer turns it off after the selected timer ends.	Restart if you want another session.
	Battery has run out.	Charge the device.
Feel too powerful, uncomfortable, or tight while massaging.	Intensity is set too high.	1. Choose lower intensity. 2. Press and hold the power button to turn off. 3. Zip down the wrap.
The controller gets hot.	It is normal for the controller to slightly heat up if the device has been used for a long time.	It is strongly advised not to use this massager for more than 60 minutes per session to avoid OVER USAGE. After turning it off, allow the controller to cool down before using again.

WARRANTY POLICY

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not void the product warranty.

CUSTOMER SERVICE

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

 Tel: 1800 MY RENPHO (1800 69 7367) (AU ONLY)
Monday-Friday 9:00AM-6:00PM (AEST)

 Email: support-au@renpho.com (AU)

 Tel: +44(0)7434-666088 (UK ONLY)
Monday-Friday 9:00AM-5:00PM (GMT)

 Email: support-eu@renpho.com (UK&EU)

*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste.

Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

1. The symbol of the crossed-out garbage can on batteries or accumulators indicates that they must not be disposed of in household waste at the end of their life.

2. You are legally obliged to return old batteries and accumulators after use. You can do this free of charge at a collection point in your area. Addresses of suitable collection points can be obtained from your city or local government.

3. Batteries may contain substances that are harmful to the environment and human health. The separate collection and recycling of old batteries and accumulators is intended to avoid negative effects on the environment and human health.

4. Please avoid the generation of waste from old batteries as far as possible, e.g. Please avoid littering public spaces by not carelessly leaving batteries or electrical and electronic equipment containing batteries.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/30/EU, Directive 2011/65/EU, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address: <https://renpho.eu/pages/compliance-data>

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