
RENPHO

User Manual

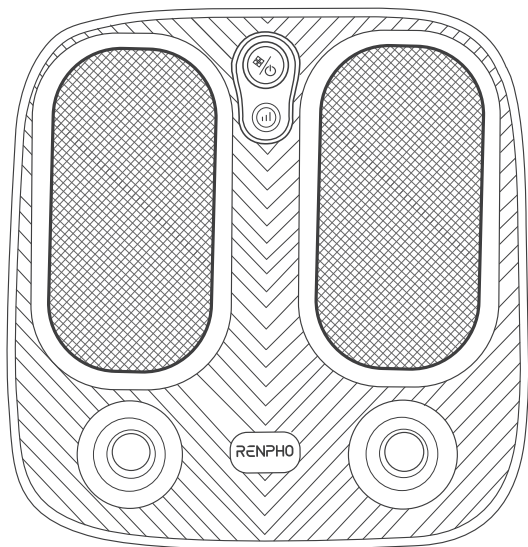TEREA SOLE

Model: R-WMF10

CONTENTS

• ENGLISH	01 ~ 09
-----------	---------

• DEUTSCH	10 ~ 18
-----------	---------

• FRANÇAIS	19 ~ 27
------------	---------

• ITALIANO	28 ~ 36
------------	---------

• ESPAÑOL	37 ~ 45
-----------	---------

IMPORTANT SAFETY INSTRUCTIONS

READ ALL INSTRUCTIONS BEFORE USING.

1. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the appliance.
3. **WARNING** – To reduce the risk of burns, fire, electric shock, or injury to persons.
 - a) This appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
 - b) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
 - c) Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
 - d) Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
 - e) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
 - f) Do not carry this appliance by supply cord or use cord as a handle.
 - g) Keep the cord away from heated surfaces.
 - h) Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
 - i) Never drop or insert any object into any opening.
 - j) Do not use outdoors.
 - k) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
 - l) To disconnect, turn all controls to the off position, then remove plug from outlet.
4. **DANGER:** Do not use while bathing or in a shower.
5. **DANGER:** Do not place or store appliance where it can fall or be pulled into a tub or sink.
6. **DANGER:** Do not place in or drop into water or other liquid.
7. **WARNING:** Never operate on a soft surface such as a bed or couch where the air openings may be blocked.

8. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
9. **WARNING:** Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
10. **WARNING:** Do not stand on or in appliance. Use only while seated.
11. The device is intended for household use only and not intended for medical use. It should not be used in treatment or alleviation of disease or abnormal physical state.

About the Battery

1. Do not expose the device or battery to excessive temperatures.
2. Be aware of the risk of terminals of the battery-operated device or battery being short- circuited by metal objects.
3. This device contains batteries that are non-replaceable. When the battery is at end of life, the device shall be properly disposed of.
4. Do not dismantle, open or shred battery.
5. Do not subject cells or batteries to mechanical shock.
6. In the event of a cell leaking, do not allow the liquid to come in contact with the skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical advice.
7. Secondary cells and batteries need to be charged before use. Always use the correct charger and refer to the manufacturer's instructions or product manual for proper charging instructions.
8. Do not leave a battery on prolonged charge when not in use.
9. After extended periods of storage, it may be necessary to charge and discharge the cells or batteries several times to obtain maximum performance.
10. **CAUTION** - The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not recharge, disassemble, heat above 100°C (212°F).

SAVE THESE INSTRUCTIONS

Battery Producer Labelling

Manufacturer: XIAMEN INVCY ELECTRONICS CO., LTD.

Address: 1F/3F-5F, BLDG 5, No.943-5, Tonglong 2nd Rd,

Torch Hi-tech Industrial Zone (Xiang'an)Xiamen 361000 Fujian P.R. China

Email: 2162923595@qq.com

ABOUT THE DEVICE

1. What's in the Box

Terea Sole × 1

USB-C Charging Cable × 1

2. Technical Specifications

Product Name: Foot Massager

Model: R-WMF10

Input: 5V  2A

Battery Capacity: 2500mAh 18.5Wh


Charging Time: About 4 h

Operation Time: > 1 h

Product Weight: 4.2 lbs / 1.9 kg

Product Size: 12 × 12.5 × 3 in / 30.4 × 31.8 × 7.5 cm

3. Overview

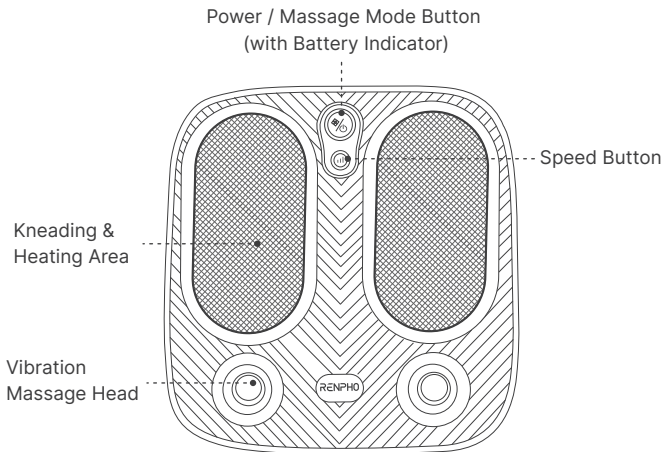USING THE DEVICE

- * Do Not Use The Device For More Than 30 Mins Each Time.
- * Do Not Stand Up While Massaging.
- * Fully Charge The Device Before Initial Use, And Do Not Use It While Charging.

Button Instruction

Button	Operation
	Press and hold to turn ON /OFF the device. Press to switch 3 massage modes: • Mode 1: Kneading + Vibration • Mode 2: Kneading + Vibration + Heat + Light • Mode 3: Intermittent Kneading & Vibration Note: Kneading and heat (30°C / 86°F) for forefeet and arches, and vibration for heels.
	Press to cycle through 3 speed levels, with the corresponding light on: • Level 1: Blue • Level 2: Orange • Level 3: Red

Note:

Stall Protection: If the device stalls for 5 seconds, it will stop operation for 3 seconds before auto-restarting.

Battery Indicator

Status	Indicator
Normal operation	Solid blue
Low battery / charging	Flashing blue
Fully charged	Solid blue

MAINTENANCE AND CLEANING

1. It is prohibited for non-professional personnel to repair the device or replace its parts.
2. Place and store the device in a cool and dry environment. Do not expose the device to direct sunlight. Keep away from humidity, heat, and/or corrosive gas.
3. Do not place any heavy objects on the device.
4. Keep the device clean. Unplug the device before cleaning.
5. Use a soft, damp cloth to wipe the main body for daily cleaning. To wipe out the stains caused by oil or cosmetics, use a soft, damp cloth with neutral detergent or special detergent for leather. Never immerse the device in water or other liquids.

TROUBLESHOOTING GUIDE

Problems	Causes	Solutions
Fail to turn on the device.	Batteries are weak or depleted.	Charge the device and restart it.
	Power Button was pressed for only a short duration.	Press and hold the Power Button to turn on the device.
Slow charging / no charging.	Use other types of charging cables.	Use the provided USB-C charging cable only.
The device does not work.	Applying excessive pressure could cause nodes to stall.	Ensure you are seated and avoid applying excessive pressure to the nodes when in use.
Insufficient heat or no heat.	The device takes time to heat up gradually to avoid discomfort caused by rapid heat.	Wait several minutes after turning on the device.
	Thick socks are worn during massage, reducing sensation.	Wear thinner socks to improve foot sensation and the overall effect of the massage.

WARRANTY POLICY

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not diminish the product warranty.

CUSTOMER SERVICE

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

TEL: +1(844) 417 0149 (US&CA)

Monday-Friday 9:00AM-4:30PM

1800 MY RENPHO (1800 69 7367) (AU ONLY)

Monday-Friday 9:00AM-6:00PM (AEST)

+44(0)7434 666088 (UK ONLY)

Monday-Friday 9:00AM-5:00PM (GMT)

Email: support@renpho.com (US&CA)

support-au@renpho.com (AU)

support-eu@renpho.com (UK&EU)

* For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

1. The symbol of the crossed-out garbage can on batteries or accumulators indicates that they must not be disposed of in household waste at the end of their life.
2. You are legally obliged to return old batteries and accumulators after use. You can do this free of charge at a collection point in your area. Addresses of suitable collection points can be obtained from your city or local government.
3. Batteries may contain substances that are harmful to the environment and human health. The separate collection and recycling of old batteries and accumulators is intended to avoid negative effects on the environment and human health.
4. Please avoid the generation of waste from old batteries as far as possible, e.g. Please avoid littering public spaces by not carelessly leaving batteries or electrical and electronic equipment containing batteries.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/30/EU, Directive 2011/65/EU, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address:
<https://renpho.eu/pages/compliance-data>

FCC Regulatory Compliance

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

ISED Regulatory Compliance

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- (1) This device may not cause interference.
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

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